

Movilo®

Therapeutic Comfort Shoe

NEW



MOVILLO

DARCO



Movilo®

Perfect for people who spend a lot of their time on their feet and value healthy, pain-free walking.

Movilo® Fly



Colours: Black, Blue, Lilac, Petrol

Movilo® Street




Colours: White, Olive

Indications

For the treatment of foot complaints:

- > Moderate foot / toe deformities such as hallux valgus, hammer and claw toes
- > Hallux rigidus
- > Metatarsalgia
- > Morton's neuroma
- > Rheumatoid arthritis
- > Polyarthritis, arthrosis, gout
- > Fasciitis plantaris / heel spur
- > Callus formation, panaritium prophylaxis
- > Diabetes mellitus without PNP / PAOD (risk group 0 according to the German Diabetes Society)

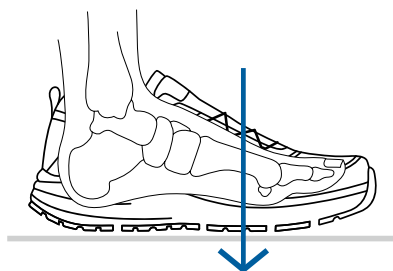

The back-shifted ball rocker facilitates the rolling off movement and reduces pressure in the forefoot area.

Features



OUTSOLE

- > Back-shifted ball rocker
 - Forefoot off-loading
- > Large abrasion-resistant rubber coating (3 mm)
 - Non-slip
- > Shock-absorbing, force-absorbing sole construction
- > Continuous, rigid shoe base reinforcement
 - Forefoot off-loading
- > Long and secure footprint
- > Lower heel-to-toe drop (11 mm)



SHAFT

- > Air-permeable
- > Reinforced toe
- > Anatomical ankle cutout
 - Ideal rearfoot fit
- > Optimally adapted heel cap
 - Good, safe heel support
- > Loop on the heel
 - Easy to put on

Model Fly

- > Loop on the shoe tongue
 - Easy to put on
- > Without liner
 - Slim
 - Soft
 - Breathable



Model Street

- > Stable upper material
- > Soft liner
 - Good padding
 - Protection against pressure and external influences
 - High wearing comfort




INSOLE SYSTEM

- > 2-part footbed with 8 mm total height
 - Anatomically moulded footbed (6 mm)
 - Moderate heel guide
 - Slight lengthwise and crosswise arch support
 - PU sole with smooth microfibre cover for good cushioning + easy foot insertion
 - Placeholder sole (2 mm)
 - For volume compensation
 - Additional shock absorption
- > Replaceable
- > Sufficient space (8 mm) for customized insoles
- > Individually adjustable



LAST "PINE"

- > Foot-friendly last shape
- > Ideal ball width
- > Toe box height 21.5 – 27 mm (depending on size)
- > Developed in cooperation with Wiest GmbH

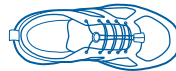

LACING

- > Elasticated laces with quick-release fastener
 - Easy to put on and take off
- > Non-elastic, flat laces (enclosed)
 - Better support
- > Additional lace hole
 - Firm fit



Article List, Sizes and Measurements




Movilo® Fly

Size	Insole in cm		Foot length in cm F	Ball girth in cm G	Toe box height in cm H	Art. No.			
	A	B				Black	Blue	Lilac	Petrol
36	24,0	9,0	22,5 - 23,1	22,6	2,15	MOFBAP36	MOFBUP36	MOFLCP36	MOFPTP36
37	24,6	9,1	23,1 - 23,7	23,0	2,20	MOFBAP37	MOFBUP37	MOFLCP37	MOFPTP37
38	25,3	9,2	23,8 - 24,4	23,4	2,25	MOFBAP38	MOFBUP38	MOFLCP38	MOFPTP38
39	26,0	9,3	24,5 - 25,1	23,8	2,30	MOFBAP39	MOFBUP39	MOFLCP39	MOFPTP39
40	26,6	9,4	25,1 - 25,7	24,2	2,35	MOFBAP40	MOFBUP40	MOFLCP40	MOFPTP40
41	27,3	9,6	25,8 - 26,4	24,6	2,40	MOFBAP41	MOFBUP41	MOFLCP41	MOFPTP41
42	28,0	9,8	26,5 - 27,1	25,0	2,45	MOFBAP42	MOFBUP42	MOFLCP42	MOFPTP42
43	28,6	10,1	27,1 - 27,7	25,4	2,50	MOFBAP43	MOFBUP43	MOFLCP43	MOFPTP43
44	29,3	10,2	27,8 - 28,4	25,8	2,55	MOFBAP44	MOFBUP44	-	MOFPTP44
45	30,0	10,3	28,5 - 29,1	26,2	2,60	MOFBAP45	MOFBUP45	-	MOFPTP45
46	30,6	10,4	29,1 - 29,7	26,6	2,65	MOFBAP46	MOFBUP46	-	MOFPTP46
47	31,3	10,5	29,8 - 30,4	27,0	2,70	MOFBAP47	MOFBUP47	-	MOFPTP47



Movilo® Street

Size	Insole in cm		Foot length in cm F	Ball girth in cm G	Toe box height in cm H	Art. No.	
	A	B				White	Olive
36	24,0	9,0	22,5 - 23,1	22,6	2,15	MOSWIP36	MOSOVP36
37	24,6	9,1	23,1 - 23,7	23,0	2,20	MOSWIP37	MOSOVP37
38	25,3	9,2	23,8 - 24,4	23,4	2,25	MOSWIP38	MOSOVP38
39	26,0	9,3	24,5 - 25,1	23,8	2,30	MOSWIP39	MOSOVP39
40	26,6	9,4	25,1 - 25,7	24,2	2,35	MOSWIP40	MOSOVP40
41	27,3	9,6	25,8 - 26,4	24,6	2,40	MOSWIP41	MOSOVP41
42	28,0	9,8	26,5 - 27,1	25,0	2,45	MOSWIP42	MOSOVP42
43	28,6	10,1	27,1 - 27,7	25,4	2,50	MOSWIP43	MOSOVP43
44	29,3	10,2	27,8 - 28,4	25,8	2,55	MOSWIP44	MOSOVP44
45	30,0	10,3	28,5 - 29,1	26,2	2,60	MOSWIP45	MOSOVP45
46	30,6	10,4	29,1 - 29,7	26,6	2,65	MOSWIP46	MOSOVP46
47	31,3	10,5	29,8 - 30,4	27,0	2,70	MOSWIP47	MOSOVP47

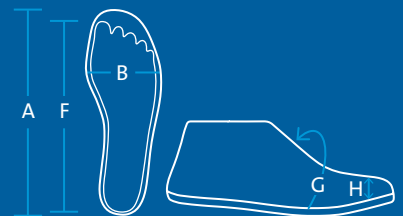

Let's go!



Find the right size for you