



MOViLO

Therapeutic Comfort Shoe

Let's go!

DARCO



Where's the rub?

- > Do you spend a lot of time on your feet and attach great importance to healthy, pain-free walking?
- > Do you have problems with your feet and find it difficult to find the right shoe?
- > Are you looking for comfort shoes with sufficient space for your individual insoles?
- > Do you want comfortable shoes that also look great?

Well then – Let's go!

Our Movilo® Fly and Movilo® Street therapeutic comfort shoes offer you plenty of room, are extremely comfortable to wear, and also look attractive. The Movilo® shoes from DARCO support you at work and in active everyday life.





Find the right
size for you



Movilo® Fly



Colours	Size
Black	36 – 47
Blue	36 – 47
Lilac	36 – 43
Petrol	36 – 47

Movilo® Street



Colours	Size
White	36 – 47
Olive	36 – 47

Legal notice

DARCO neither practices in the medical field, nor does it make recommendations for specific patient use. The specialist medical staff is responsible in each individual patient case for determining the appropriate devices and for carrying out the suitable treatment. DARCO accepts no responsibility for this.



Features



OUTSOLE

- > Back-shifted ball rocker
 - Forefoot off-loading
- > Large abrasion-resistant rubber coating (3 mm)
 - Non-slip
- > Shock-absorbing, force-absorbing sole construction
- > Continuous, rigid shoe base reinforcement
 - Forefoot off-loading
- > Long and secure footprint
- > Lower heel-to-toe drop (11 mm)



INSOLE SYSTEM

- > 2-part footbed with 8 mm total height
 - Anatomically moulded footbed (6 mm)
 - Moderate heel guide
 - Slight lengthwise and crosswise arch support
 - PU sole with smooth microfibre cover for good cushioning + easy foot insertion
 - Placeholder sole (2 mm)
 - For volume compensation
 - Additional shock absorption
- > Replaceable
- > Sufficient space (8 mm) for customised insoles
- > Individually adjustable



SHOE LAST "PINE"

- > Foot-friendly last shape
- > Ideal ball width
- > Toe box height 21.5 – 27 mm (depending on size)
- > Developed in cooperation with Wiest GmbH



SHAFT

- > Air-permeable
- > Reinforced toe
- > Anatomical ankle cutout
 - Ideal rearfoot fit
- > Optimally adapted heel cap
 - Good, safe heel support
- > Loop on the heel
 - Easy to put on


Model Fly

- > Loop on the shoe tongue
 - Easy to put on
- > Without liner
 - Slim
 - Soft
 - Breathable



Model Street

- > Stable upper material
- > Soft liner
 - Good padding
 - Protection against pressure and external influences
 - High wearing comfort

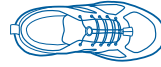

LACING

- > Elasticated laces with quick-release fastener
 - Easy to put on and take off
- > Non-elastic, flat laces (enclosed)
 - Better support
- > Additional lace hole
 - Firm fit



Indications

For the treatment of foot complaints:

- > Moderate foot / toe deformities such as hallux valgus, hammer and claw toes
- > Hallux rigidus
- > Metatarsalgia
- > Morton's neuroma
- > Rheumatoid arthritis
- > Polyarthritis, arthrosis, gout
- > Fasciitis plantaris and heel spur
- > Callus formation, panaritium prophylaxis
- > Diabetes mellitus without PNP / PAOD (risk group 0 according to the German Diabetes Society)



DARCO (Europe) GmbH | Gewerbegebiet 18 | 82399 Raisting | Germany
Tel. +49 8807 9228-0 | Fax -22 | info@darco.de | www.darco.de/en



2024-07-18