

Movilo®

Therapeutic Comfort Shoe

NEW



MOVILLO

DARCO



Movilo®

Perfect for people who spend a lot of their time on their feet and value healthy, pain-free walking.

Movilo® Fly



Colours: Black, Blue, Lilac, Petrol

Movilo® Street




Colours: White, Olive

Indications

For the treatment of foot complaints:

- > Moderate foot / toe deformities such as hallux valgus, hammer and claw toes
- > Hallux rigidus
- > Metatarsalgia
- > Morton's neuroma
- > Rheumatoid arthritis
- > Polyarthritis, arthrosis, gout
- > Plantar fasciitis and heel spur
- > Callus formation, panaritium prophylaxis
- > Diabetes mellitus without PNP / PAOD (risk group 0 according to the German Diabetes Society)

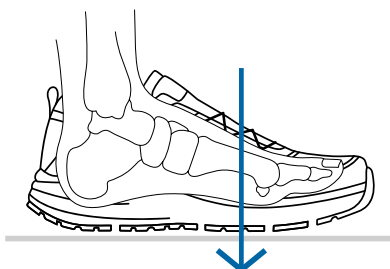

The back-shifted ball rocker facilitates the rolling off movement and reduces pressure in the forefoot area.

Features



OUTSOLE

- > Back-shifted ball rocker
 - Forefoot off-loading
- > Large abrasion-resistant rubber coating (3 mm)
 - Non-slip
- > Shock-absorbing, force-absorbing sole construction
- > Continuous and reinforced rigid shoe base
 - Forefoot off-loading
- > Long and secure footprint
- > Lower heel-to-toe drop (11 mm)



SHAFT

- > Air-permeable
- > Reinforced toe
- > Anatomical ankle cutout
 - Ideal rearfoot fit
- > Optimally adapted heel cup
 - Substantial and cushioned heel support
- > Loop on the heel
 - Easy to put on

Model Fly

- > Loop on the shoe tongue
 - Easy to put on
- > Without liner
 - Slim
 - Soft
 - Breathable



Model Street

- > Stable upper material
- > Soft liner
 - Good padding
 - Protection against pressure and external influences
 - Comfortable to wear




INSOLE SYSTEM

- > 2-part footbed with 8 mm total height
 - Anatomically moulded footbed (6 mm)
 - Moderate heel guide
 - Slight lengthwise and crosswise arch support
 - PU sole with smooth microfibre cover for good cushioning + easy foot insertion
 - Placeholder sole (2 mm)
 - For volume compensation
 - Additional shock absorption
- > Replaceable
- > Sufficient space (8 mm) for customised insoles
- > Individually adjustable



SHOE LAST MOVILLO®

- > Foot-friendly last shape
- > Ideal ball width
- > Toe box height 21.5 – 27 mm (depending on size)
- > Developed in cooperation with Wiest GmbH

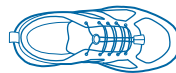

LACING

- > Elasticated laces with quick-release fastener
 - Easy to put on and take off
- > Non-elastic, flat laces (spare set enclosed)
 - Better support
- > Additional lace hole
 - Secure fit




Article List, Sizes and Measurements

Movilo® Fly

Size	Insole in cm		Foot length in cm F	Ball girth in cm G	Toe box height in cm H	Art. No.			
	A	B				Black	Blue	Lilac	Petrol
36	24.0	9.0	22.5 - 23.1	22.6	2.15	MOFBAP36	MOFBUP36	MOFLCP36	MOFPTP36
37	24.6	9.1	23.1 - 23.7	23.0	2.20	MOFBAP37	MOFBUP37	MOFLCP37	MOFPTP37
38	25.3	9.2	23.8 - 24.4	23.4	2.25	MOFBAP38	MOFBUP38	MOFLCP38	MOFPTP38
39	26.0	9.3	24.5 - 25.1	23.8	2.30	MOFBAP39	MOFBUP39	MOFLCP39	MOFPTP39
40	26.6	9.4	25.1 - 25.7	24.2	2.35	MOFBAP40	MOFBUP40	MOFLCP40	MOFPTP40
41	27.3	9.6	25.8 - 26.4	24.6	2.40	MOFBAP41	MOFBUP41	MOFLCP41	MOFPTP41
42	28.0	9.8	26.5 - 27.1	25.0	2.45	MOFBAP42	MOFBUP42	MOFLCP42	MOFPTP42
43	28.6	10.1	27.1 - 27.7	25.4	2.50	MOFBAP43	MOFBUP43	MOFLCP43	MOFPTP43
44	29.3	10.2	27.8 - 28.4	25.8	2.55	MOFBAP44	MOFBUP44	-	MOFPTP44
45	30.0	10.3	28.5 - 29.1	26.2	2.60	MOFBAP45	MOFBUP45	-	MOFPTP45
46	30.6	10.4	29.1 - 29.7	26.6	2.65	MOFBAP46	MOFBUP46	-	MOFPTP46
47	31.3	10.5	29.8 - 30.4	27.0	2.70	MOFBAP47	MOFBUP47	-	MOFPTP47

Movilo® Street

Size	Insole in cm		Foot length in cm F	Ball girth in cm G	Toe box height in cm H	Art. No.	
	A	B				White	Olive
36	24.0	9.0	22.5 - 23.1	22.6	2.15	MOSWIP36	MOSOV36
37	24.6	9.1	23.1 - 23.7	23.0	2.20	MOSWIP37	MOSOV37
38	25.3	9.2	23.8 - 24.4	23.4	2.25	MOSWIP38	MOSOV38
39	26.0	9.3	24.5 - 25.1	23.8	2.30	MOSWIP39	MOSOV39
40	26.6	9.4	25.1 - 25.7	24.2	2.35	MOSWIP40	MOSOV40
41	27.3	9.6	25.8 - 26.4	24.6	2.40	MOSWIP41	MOSOV41
42	28.0	9.8	26.5 - 27.1	25.0	2.45	MOSWIP42	MOSOV42
43	28.6	10.1	27.1 - 27.7	25.4	2.50	MOSWIP43	MOSOV43
44	29.3	10.2	27.8 - 28.4	25.8	2.55	MOSWIP44	MOSOV44
45	30.0	10.3	28.5 - 29.1	26.2	2.60	MOSWIP45	MOSOV45
46	30.6	10.4	29.1 - 29.7	26.6	2.65	MOSWIP46	MOSOV46
47	31.3	10.5	29.8 - 30.4	27.0	2.70	MOSWIP47	MOSOV47

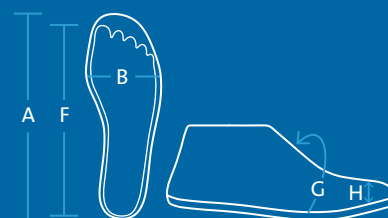

Find the right size for you